

T-Connect Quick Start Guide



1


Take the T-Top machine out of the box

2



Turn on The T-Top machine

3

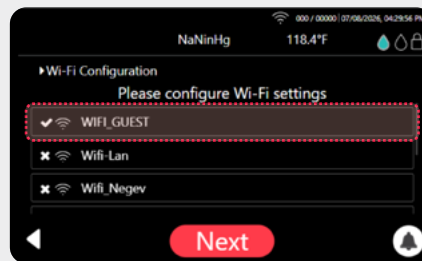


In the Startup Wizard, connect to Wi-Fi during the wizard process

4




Upon completion press the 3 lines in the top left corner of the screen

4

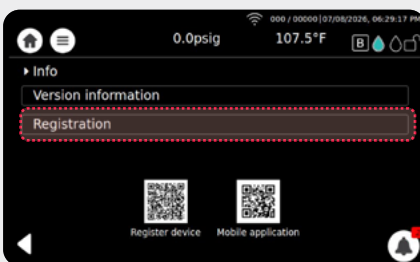

Select Registration

5

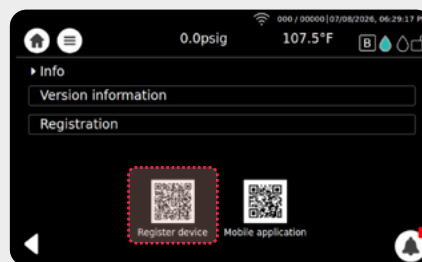

Scan QR Code to register your device

6



An e-mail will be sent to you to log into T-Connect software, the name of your clinic and machine will appear upon log in



For immediate service or technical support contact our support team
1-800-930-0452 or techsupport@tuttnauerusa.com