TIVA2 Washers

Thermal High Disinfector Line

Dental and medical practices need strict cleaning, disinfection, and sterilization of instruments before they can be used on patients. Thermal High Disinfection employs moist heat above 194° for a minimum of 5 minutes to kill bacteria and viruses in the washing phase. The TIVA2 line provides:

- Thermal High Disinfection to achieve a complete decontamination cycle.
- Advanced control system and an easy to use interface

TIVA2-H TD

Under Counter Washer

TIVA2 TD

Counter Top Washer





High end washers for the decontamination of complex instruments

TIVA2 Thermal High Disinfector washers provide a flexible turnkey solution for washing, disinfecting and drying diverse medical instruments for a complete decontamination cycle.

a minimum of 5

minutes

phase



alkaline detergents



TIVA2/TIVA2-H TD are fully-automatic thermal high disinfectors for all practices.

Power switch

Control panel

Air suction filter

Instrument rack

Washing chamber door

Detergent containers

Drawer door with lock

The front and back panel of the TIVA2 and TIVA2-H are identical

Features:

- Touch panel, color graphic display and real-time monitoring
- A total of 40 programs: Preset and programmable cycles
- Air drying prefilter + HEPA H14 filter
- Integrated water softener
- 2 dosing pumps with flowmeter for alkaline detergent and neutralizer (additional pump optional for TIVA2-H TD)
- USB/Printer for data documentation
- Provides complete & total evacuation of previous water (ETS)

Utility Requirements:

Water Connections

- Cold water connection for TIVA2 (default)
- Connection for demineralized water
- Integrated Water Softener with Salt Box
- R.O. water system is required
- Failure to incorporate an RO system with our TIVA unit will void your warranty

Electrical Connections


- 208+240V 1 phase 60Hz

Dimensions: TIVA2-H TD Under Counter Washer

- (WxDxH) 23.4"x21.3"x32.2"
- Usable Chamber Volume: 65L

TIVA2 TD Counter Top Washer

- (WxDxH) 23.4"X21.3"X23.6"

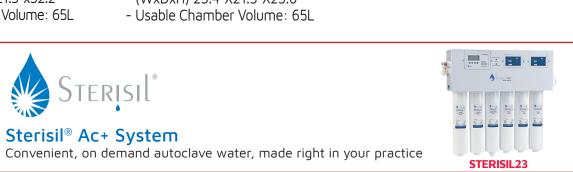

NOTE:

T-Clean Tuttnauer Line of Detergents & Neutralizers -User friendly for use in your Washer

T-Clean TIVA Detergent



TD-1L



T-Clean TIVA Neutralizer:



T-Clean Prep Foam PF0006 500 ml with sprayer









25 Power Drive, Hauppauge, NY 11788

Tel: +800 624 5836, +631 737 4850, Fax: +631 737 0720 E-mail: info@tuttnauerUSA.com, www.tuttnauerUSA.com



Visit our blog: tuttnauerusa.com/blog







Back connection panel

